BBWC TURBO POWER

(type 921500)

PAGE

WARNING / LET OP / WARNUNG / ATTENTION / ATTENZIONE / AVISO /	2
Gebruiksaanwijzing gezichtsbruiner	4
User's manual facial tanner	9
Bedienungsanleitung Gesichtsbräuner	14
Mode d'emploi du bronzeur visage	19
Istruzioni per l'uso della lampada abbronzante per il viso	24
Manual de Instrucciones de la lámpara de bronceado facial	29
Zijpanelen, sidepanels, Seitenteile, panneaux lateraux Panelli laterali	35
Declaration of Conformity	36

Rev.01.06.2024

**WARNING / LET OP / WARNUNG / ATTENTION /
ATTENZIONE / AVISO /**

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GB : This is a very powerfull professional product. We advise you to consult your doctor prior to use. Please keep the minimum distance and do not exceed the recommended exposure times as described in the manual.

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WARNING:

It is possible that the lamps don't make any contact because they moved a little during transportation. This problem can be solved by turning the lamps a little. Beware: the lamps are always switched on in combination with another lamp. So there are always 2 lamps, or a plural of 2 not working. Does drawing the one not work try the other!

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Contact

BBWC
Kapel Avezaath
THE NETHERLANDS
www.bbwc.nl

User's manual facial tanner

Please read these instructions carefully before assembling and using the sun lamp GB9215. Make sure you only use this appliance in accordance with the recommended treatment times, tanning intervals, and tanning distance.

Power supply


This sun lamp should only be connected to a 230V 50Hz socket. Never place the sun lamp in a bathroom, sauna, or any other damp room.

Protective goggles

Ultraviolet radiation is harmful to the eyes. A pair of protective goggles, which must be worn with every use, has therefore been included with this sun lamp. Never use "regular" sunglasses instead of the protective goggles. Even shutting the eyes provides inadequate protection. A spare pair is also included.

Operation

Start the facial tanner by opening the two side panels by pushing button **A** (see drawing **1**). Panels will open themselves by a built in spring. Open the panels till they cannot go any further **1 B**. DON'T FORCE SINCE THAT MAY CAUSE DAMAGE. The panels should stay in a certain angle to give maximum side tanning on the face. Turn the facial in the right position towards the face. The chosen position can be fixed by turning the knob on the right side of the facial (turn to the left **1 E**). Plug the mains lead into the socket and set the timer **1 C**

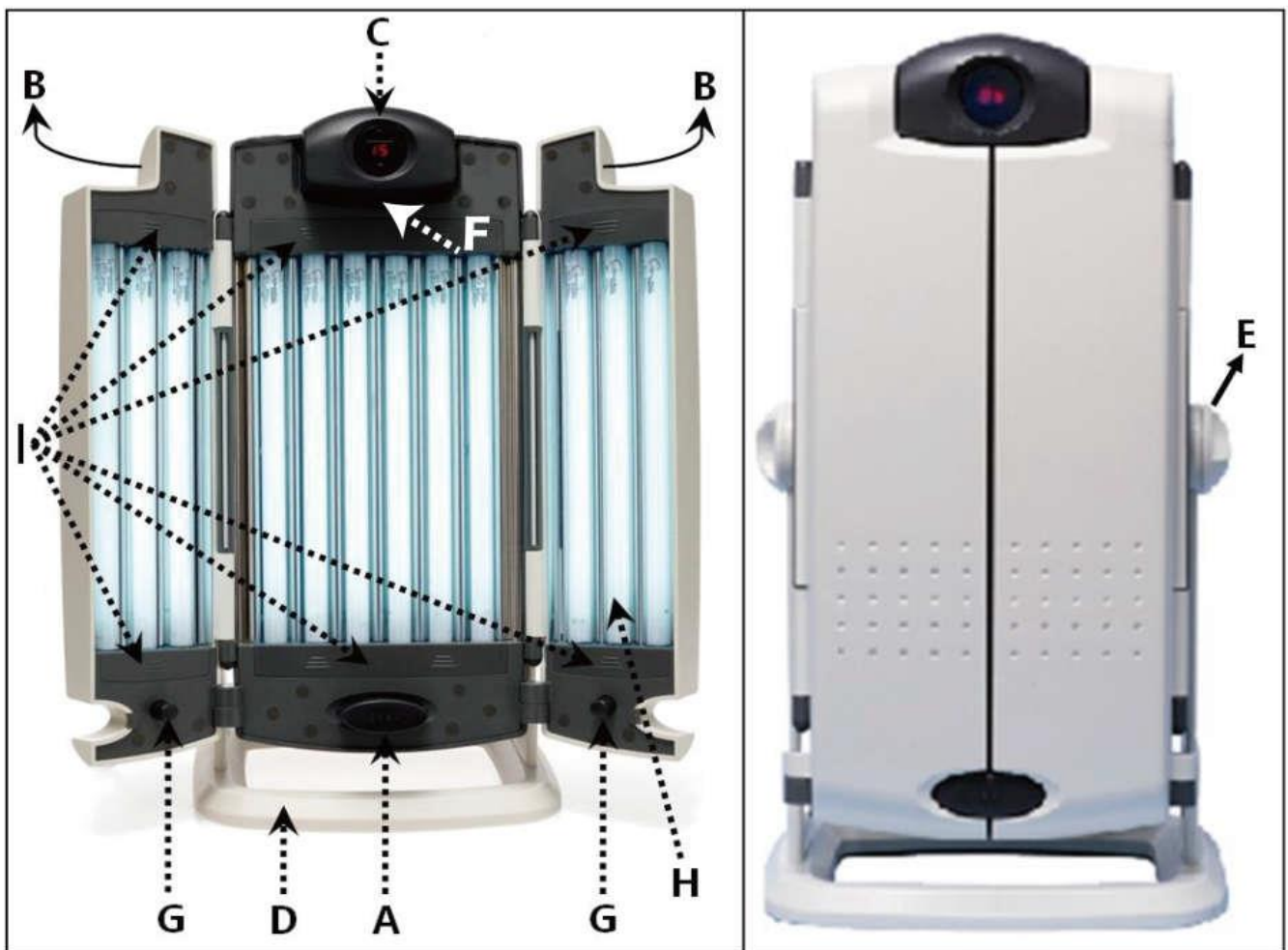DRAWING 1

Timer

Plug the main lead into the socket. The timer display will now show 00 (blinking). Set the timer to the desired time by pressing the button above the display. Each press on the button means 5 minutes. You can choose time from 5 to 30 minutes. After you have selected the time and you do not touch any other button anymore the lamps will start burning after 3 seconds. Once the preset time has run out, the facial tanner will switch off automatically. If you wish to stop using the facial tanner before the time simply press the stop button below the display. During operation the timer will give a bleep every minute. The time left will be shown on the display.

Maintenance, cleaning and technical information

First remove the plug from the socket! Wipe the outside of the sun lamp with a slightly damp cloth. Do not use an abrasive cloth. Make sure that no water gets into the sun lamp. The light tubes can be dusted with a dry cloth or a soft brush. The reflector plate should only be dusted with a very soft brush (if necessary) from time to time. The lamps have a very long operating life of about 300 hours. This facial tanner has been built according to all European safety standards, which allows it to carry the CE-mark. Insulation class II, UV radiation type 4, mains voltage 230 Volts, 50 Hz, power 200 Watts.

General information about tanning

You can tan safely all year round underneath this sun lamp if you follow these guidelines correctly. Our skin is tanned by invisible ultraviolet radiation from the sun. This radiation darkens the pigment already present in our skin and stimulates the formation of more pigment known as melanin. The production of new melanin cells and the darkening of these make the skin stay tanned longer. There are a number of different types of ultraviolet radiation. Ultraviolet A radiation (UVA), which colors the pigment cells brown, and ultraviolet B radiation (UVB), which stimulates the production of melanin cells, but also causes sunburn, are the most well-known types. How quickly and how deep you tan depends on your skin type. The skin type also determines how long one can (safely) stay in the sun. This is true for both the real sun as well as the sun lamp. In order to make the best use of your sun lamp, first check your skin type using the following table. Further on in these instructions, you will find the different skin types in the tables with the recommended tanning times. If your skin falls into the type 1 or 2 category, we recommend you do not use the facial tanner at all.

SKINTYPE

Characteristics	Type 1	Type 2	Type 3	Type 4	Type 5	Type 6
Eye color?	Light blue, grey or green	Blue, grey or green	Dark blue, grey or green	Dark mostly brown	Dark brown	Dark brown
Hair color?	Red, light blond	blond	Dark blonde, brown	Dark brown	Dark brown, black	Black
Skin color?	Very pale	pale	Slightly tinted	Tinted	Dark Tinted	Very dark
Freckles	Lots very often	sometimes	Rarely	None	None	None
Tanning experience	No tan at all	Slowly difficult	Good	Very good and deep	Extremely good	Almost no effect
Sunburns	Always	Easily	Occasionally	Rarely or never	Only by exception	Never

The tanning program for facial Eurosolar GB9215

Before you begin, read all the instructions!

exposure	Exposure schedule (in minuts) maximum				
	skintype 2	skintype 3	skintype 4	skintype 5	skintype 6
1	5	5	5	5	5
2	5	10	10	15	15
3	5	10	10	15	15
4	5	10	15	25	25
5	5	10	15	25	30
6	5	15	20	30	35
7	5	15	20	30	35
8	5	15	20	30	35
9	5	20	25	35	35
10	5	20	25	35	35

Between two exposures there should be at least 48 hours! After five exposures there should then be an interval of at least three days, then commence five exposures once again with a minimum 48 hours between exposures

Maximum number of exposures a year : 40

These recommended hours only apply for use with the recommended light tubes. Determine your skin type using the table. The maximum number of minutes for one tanning session can be found there.

Keep a minimum of 25 cm distance between your face and the lamp. A greater distance is allowed although the tanning effect will diminish drastically. Always use the protective goggles! Should your skin become red and/or tight and sensitive, it is recommended you decrease the treatment time by 25%. After a treatment of 10 exposure days, a maintenance tanning time of 30 minutes a week is sufficient to keep your tan but don't exceed the maximum number of treatments. Unplug the lead after each use. Let the lamps cool down for a while. If you wish you can close the panels by pushing them slightly. The panels will be locked automatically. **BEFORE CLOSING THE PANELS LAMPS SHOULD BE TURNED OFF!** If this by accident is forgotten the two safety switches will automatically shut down the operation. However don't use this as a regular procedure!

WARNINGS!!

- Read instructions carefully.
- UV radiation may cause damages of the eyes and skin like skin ageing and possibly skin cancer. Certain medicine and cosmetics may increase the sensitivity.
- Persons not using the solarium, especially children, should not be present in the close proximity when the unit is in operation.

UV applicances must not be used by:

- persons under the age of 18 years;
- Persons who tent to freckle;

- Persons with the natural pale or reddish skin; Persons suffering from sunburn; Persons with abnormal skin discolorations or disorders;
- Persons having more than 16 moles (2 mm or larger in diameter) on the body;
- Persons having atypical moles (atypical moles are considered as asymmetrical moles with a diameter more than 5 mm with different pigmentation and abnormal lines; in case of doubt, seek medical advice);
- Persons not able to tan at all or not able to tan without burning when exposed to the sun;
- Persons who burn easily when exposed to the sun;
- Persons who repeatedly experienced severe sun burns during their childhood;
- Persons suffering from or previously suffering from skin cancer or predisposed to skin cancer
- Persons with the next of kin having skin cancer.
- Persons with disorders arising out of increased sensitivity to light and undergoing treatment of such disorders; Persons with administered medication which increases or may increase sensitivity to light. Exposures should not exceed a person's minimal erythema dose (MED). If erythema (skin reddening) appears hours after any exposure, further exposure should cease. After one week, exposures could be restarted from the beginning of the schedule of exposure.
- In case that after the first tanning session, secondary effects like itching should occur within 48 hours, please consult your physician before the next tanning session. For the general advice concerning the recommended distance to the UV light source please refer to page 10.
- For the general advice concerning the recommended tanning times and schedule please refer to page 10.
- This device has not been designed for being used by persons (including children) that are physically or mentally handicapped or have sensory disabilities. Moreover, it is not allowed to use this device without sufficient experience and/or knowledge, unless these persons will be supervised by persons responsible for their security or in case they have been instructed how to use this device.
- Children are not aware of the dangers which may occur when operating electrical appliances. Children must therefore never be allowed to operate electrical appliances without supervision.
- Before you plug the power cord into the socket, please take care that the voltage on the label is compatible with your main power supply.
- Do not use the appliance in bath rooms, damp rooms or near swimming pools.
- Do not switch on the unit in any area where there may be explosive vapors.
- The device has to be connected to the mains supply according to the technical data sheet.
- The air-intake slots and hot-air exhaust openings on the cooling fans must not be blocked or covered, as the system may otherwise overheat.
- Special materials, which fade in case of light (for example painting), should also not be too long in the UV-light of the appliance.
- It is possible that lamps will flicker in new installations, or that moving shadow effects occur. These have no influence on the tanning performance.


CHANGING THE LAMPS

If tanning results are decreasing a lot (mostly after 300 hours use) you may change the lamps. USE **COSMOLUX S bleu 15W/10335** and **S pink 15W/17035** or an alternative 15W lamp with UV code range 15-0-(18-24) / (7,7-10,2)

1. make sure the facial tanner is switched off and the power cord has been taken out of the socket
2. remove the covers (**I** on drawing **1**) by pushing them in the direction indicated on the plastic
3. remove the tubes by turning them slightly and carefully
4. put in the new tubes and fix them by turning them into opposite direction
5. put the covers back in place

Parts and functions see drawing 1

- A Release button to open the panels
- C Timer with knob
- D Stand
- E Turning point with fixation knobs
- F Magnetic safety switches (turning off the lamps when the panels are accidentally closed during operation. (NO STANDARD PROCEDURE). See instruction manual.
- G Spring buttons to release the panels easier
- H Lamps 6xCosmolux S bleu 15W/10335 and 6xCosmolux S pink 15W/17035
- I Lamp covers

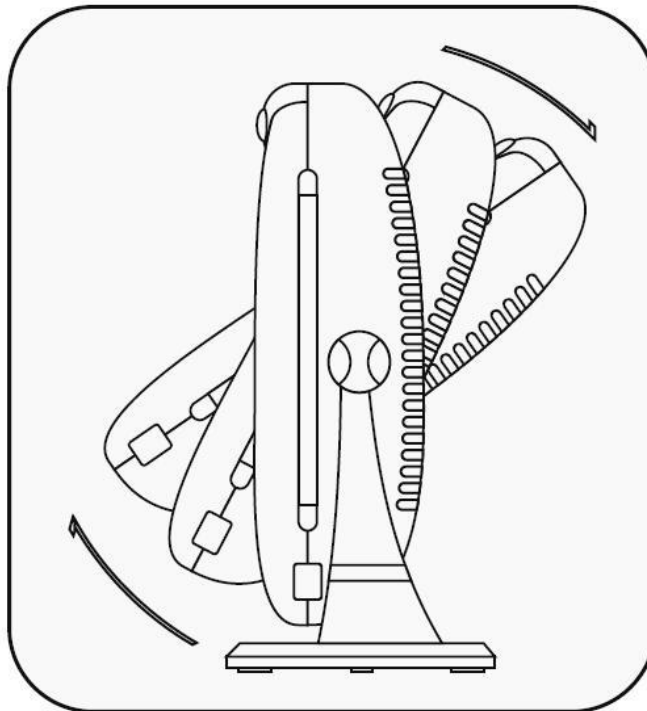Different positions to adjust with fixation knob E

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IMPORTANT INFORMATION ON THE SIDE PANELS

PLEASE BE INFORMED THAT THE SIDE PANELS HAVE NO MECHANICAL ASSISTANCE TO KEEP THEM OPEN. IF THE UNIT IS USED IN THE CORRECT POSITION AS DESCRIBED IN THE INSTRUCTIONS UNDER CHAPTER "OPERATION" ON PAGE 1. THE PANELS WILL STAY OPEN. TO ADJUST THE ANGLE OF THE UNIT USE THE RIGHT HAND KNOB TO TILT IT BACKWARDS. DO NOT ATTEMPT TO USE THE LEFT HAND KNOB AS IT IS FIXED AND UNMOVABLE

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Declaration of Conformity

The BBWC 9215 Turbo Power facial tanner meets all current and valid European Norms with regard to electrical safety and environment. BBWC BV confirms that this product complies with LVD and EMC as well as RoHs. Consequently the BBWC 9215 Turbo Power complies with all European norms below:

LVD

IEC 60335-2-27: 2019

IEC 60335-1: 2010 + Cor1:2010, Cor2:2011 + AMD1:2013, Cor1:2014+AMD2:2016, Cor1:2016

EN 60335-1: 2012 +A11:2014+A13:2017+A1:2019+A2:2019+A15:2021

EN 62233:2008

EN 60335-2-27:2013

EMC

EN IEC 55014-1: 2021

EN IEC 55014-2: 2021

EN IEC 61000-3-2: 2019+A1:2021

EN61000-3-3: 2013+A1:2019

RoHs

RoHs 2011/65/EU & 2015/863/EU

IEC 62321-3-1:2013;IEC 62321-4:2013+AMD1:2017;

IEC 62321-5:2013;IEC 62321-7-1:2015;

IEC 62321-7-2:2017;IEC 62321-6:2015;

IEC 62321-8:2017

Based on STC reports :LVD-D224656VOC 2022, EMC-D226630VOC 2022 and DTI20242513R 2024 This product carries the CE mark fully justified

BBWC

4013CN

Netherlands

www.bbwc.nl

email: info@bbwc.nl